

C a t a l o g

*Femina*  
by  EDM

[www.edm-imaging.com](http://www.edm-imaging.com)

**EDM**  
medical  
imaging

*Une autre vision  
de l'imagerie*





## For over 30 years,

EDM has been providing its customers with customized solutions that meet their requirements for accessories and consumables in Imaging. The Fémina range, dedicated to gynecology, is here combined with a wide selection of complementary products specific to your needs, whether you are a midwife or a gynecologist.

### Our dedicated team for office practices is here to listen to you:



By phone:

**+33 (0)1 30 11 83 00**

Monday to Friday from 9:00 am to 6:00 pm,  
Friday from 9:00 am to 5:00 pm



By fax:

**+33 (0)1 39 35 75 25**



On e-shop:

**[www.shop.edm-imaging.com](http://www.shop.edm-imaging.com)**



By mail:

**1, rue des Boisseliers - 95330 Domont  
France**

# SUMMARY

## DISPOSABLE INSTRUMENTATION



Functional exploration .....	p 4
Smears - Cytological sampling .....	p 7
Biopsy and histological sampling .....	p 8

## REUSABLE INSTRUMENTATION



Functional exploration .....	p 12
Stainless steel instrumentation .....	p 13
Contraception .....	p 15
Biopsy and sampling .....	p 16
Instrument disinfection .....	p 18
Sterilization .....	p 20

## WOMEN'S IMAGING



Hysterosalpingography .....	p 22
Imaging consumables .....	p 23
Ultrasound gel .....	p 25
Ultrasound probe covers .....	p 29
Ultrasound probe disinfection .....	p 30

## CONSUMABLES



Examination drape .....	p 34
Protection and attire .....	p 35
Gloves and finger cots .....	p 37
Care products .....	p 39
Surface disinfection .....	p 44
Hand hygiene .....	p 45

## EQUIPMENT



Furniture .....	p 47
Media .....	p 48
Batteries .....	p 51




## CUSCO EDM SPECULUM

*Ideal for your gynecological examinations. Wide range of sizes for every patient!*

- Enhanced visibility: speculum transparency
- Adaptation to 1 unit morphology: 6 different sizes (from 10.9 to 12 cm)
- Improved access to the cervix
- Atraumatic use: device with locking mechanism
- Avoids risks of cross-contamination

### Technical specifications:

- Class I sterile
- Aseptic crystal polystyrene
- Ethylene oxide sterilization
- Latex-free
- Sterile individual packaging

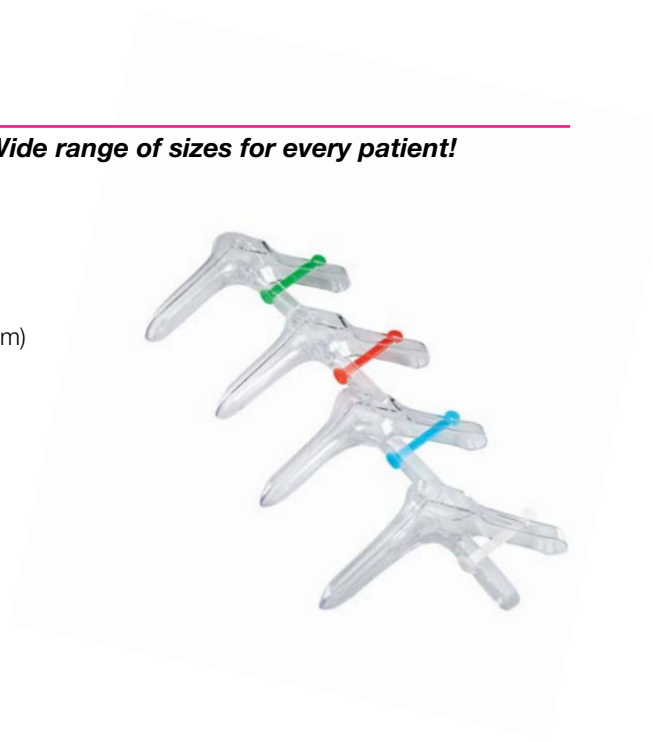

SKU	Size	Color	Diam.	Length	Qty
SPI204	XXS	●	15 mm	109 mm	100
SPI205	XS	○	20 mm	110 mm	100
SPI208	S	●	24 mm	109 mm	100
SPI206	M	●	26 mm	109 mm	100
SPI207	L	●	30 mm	112 mm	100
SPI209	XL	●	31 mm	120 mm	100

## LEGUEU TYPE FINGER

Used during medical examinations in gynecological exploration

**Technical specifications:**

- One size fits all: length 24 cm x width 20 cm
- Single use
- Made of opaque polyethylene



SKU	Description	Qty
DP422	Non-sterile - 2 fingers	100
DS426	Non-sterile - 2 fingers, individually packaged	100



## EDM LUBRICATING GEL

Exclusivity

Non-greasy, transparent, and non-irritating. Fragrance-free, alcohol-free, and hypoallergenic

- Water-based formula
- Gentle on mucous membranes
- Compatible with medical instrumentation
- Non-dripping and non-sticky
- Non-greasy
- Available in 3 packaging options



SKU	Description	Qty
KLY5	Sterile 5g	100
KLY5-S	Sterile 5g packet	100

SKU	Description	Qty
KLY42	42 g tube	6
KLY82	82 g tube	12
KLY006	250 ml dispenser	1

## Contraception

### SUBCUTANEOUS GYNECOLOGICAL IMPLANT REMOVAL KIT

Care for the incision made during the removal of the contraceptive implant

- 2 non-woven compresses 7.5 x 7.5 cm
- 1 drape 50 x 60 cm with adhesive aperture
- 1 syringe 3 ml
- 1 orange needle (25 G)
- 1 scalpel N° 15
- 1 adhesive strip 10 x 120 mm
- 1 elastic bandage 3 m x 5 cm
- 1 curved Halstead mosquito forceps 12 cm
- 1 rigid blister pack with 3 compartments



SKU	Description	Qty
EDMM50202	Subcutaneous Gynecological Implant Removal Set	5

**Technical specifications:**

- Class I sterile
- Ethylene oxide sterilization





SKU 1403.0003



## ENDOEXOCOL BRUSH

**Allows sampling of both the endocervix and the exocervix in a single rotating motion**

- All cells are collected in a single motion, allowing for simultaneous sampling of endocervical and exocervical cells
- Flexible and detachable brush, ideal for liquid-based cytology
- Thanks to its ergonomics and flexibility, the brush perfectly conforms to the contours of the cervix



INDICATIONS Sampling of the endoexocervix

### Technical specifications:

- Class I sterile
- Ethylene oxide sterilization

SKU	Description	Qty
1403.0003	Sterile endoexocervical brush	100




## SWAB

**Allows for cytological swabbing of the endocervix**



INDICATIONS Vaginal and endocervical sampling

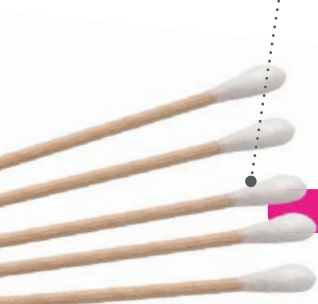

SKU 1403.0008

SKU	Description	Length	Qty
1403.0008	Swabs	150 mm	100



SKU 1403.0004

## CERVICAL BRUSH

For the collection and improvement of the endovaginal smear

**Technical specifications:**

- Gentle and effective sampling
- Its bulb reduces the risk of cervical perforation during the examination.



Endocervical smears, or any bacteriological, mycological, and/or parasitological sampling



SKU 1403.0007

## AYRE SPATULA

For sampling during uterine smears



Exocervical sampling

**Technical specifications:**

- Soft wooden spatula
- Non-sterile

SKU	Description	Length	Qty
1403.0004	Brush Ø 3.5 mm	190 mm	100
1403.0005	Brush Ø 5.5 mm (with bulb)	200 mm	100



SKU EDM211751

## CARE SWAB

For the collection and improvement of the endovaginal smear

SKU	Description	Length	Qty
EDM211751	Care swab	150 mm	100

**Technical specifications:**

- Class I - non-sterile
- 15 cm polypropylene shaft.
- 3 cm cotton tip.

## STERILE PLASTIC CUP

For the application of ointments

**Technical specifications:**

- Class I - sterile
- Ethylene oxide sterilization
- Single packaging
- Capacity: 120 ml
- Diameter: 80 mm



SKU	Description	Qty
EDM174617	Sterile plastic cup 120 ml	100

## ENDOMETRIAL ASPIRATION CURETTE

### Endometrial biopsy by aspiration



- Non-traumatic for the patient thanks to the rounded tip (SKU 1403.0010)
- Better sampling of all areas of the uterine cavity, even difficult-to-r1 unit areas, thanks to its flexibility
- More effective aspiration with the 4-hole curette
- Graduation from 2 to 12 cm



- Histological and cytological diagnosis of the endometrium
- Detection of endometrial hyperplasia, endometrial cancer, and other gynecological conditions
- Determination of response to estrogen-based hormone replacement therapy
- Bacterial or microscopic culture during uterine menstrual sample collection
- Detection of pathologies leading to infertility

#### Technical specifications:


- Z-type curette with 3 or 4 aspiration holes
- Depth marker located at 2, 4, 6, 8, 10, 12
- Class I sterile
- Ethylene oxide sterilization

SKU	Description	Total Length	Qty
1403.0009	Curette with 4 aspiration holes	275 mm + Ø 3 mm	100
1403.0010	Curette with 3 aspiration holes	285 mm + Ø 3 mm	100



## ENDOMETRIAL ASPIRATION CURETTE WITH SYRINGE

- Facilitated handling thanks to the curved and graduated curette
- 3 mm semi-rigid curette to minimize the risk of perforation
- Allows for multiple samplings during a single intervention
- Designed to improve the efficiency and reliability of sampling and reduce the patient's physical discomfort
- At the end of sampling, the in-place curette can be used as an injection route for contrast medium and for performing hysterosonography

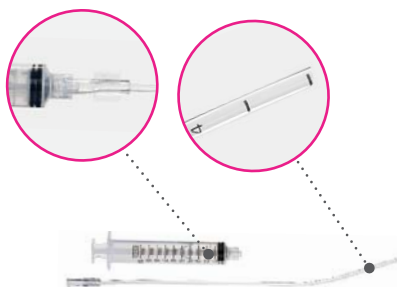

- Histological and cytological sampling of the endometrium
- Recommended for pre- and post-menopausal women for diagnostic purposes

#### Technical specifications:

- Y-shaped device consisting of a 10 cc syringe with a stopper spring to prevent reflux, and a 3 mm diameter curette
- Class I sterile
- Ethylene oxide sterilization

SKU	Volume	Total Length	Diameter	Qty
1403.0011	10 cc	235 mm	3 mm	25



Biopsy and histological sampling

**VULVAR BIOPSY PUNCH**

Ensures a core biopsy

- Multiple sizes: 2, 3, 4, 5, 6, and 8 mm
- Secure grip thanks to the ribbed handle
- Easy identification with size engraved on 1 unit punch



- Diagnosis of certain anomalies of the vulvoperineal region




SKU	Description	Diameter	Qty
2302.3000	Biopsy punch	1 mm	25
2302.3100	Biopsy punch	2 mm	25
2302.3200	Biopsy punch	3 mm	25
2302.3300	Biopsy punch	4 mm	25
2302.3400	Biopsy punch	5 mm	25
2302.3500	Biopsy punch	6 mm	25
2302.3700	Biopsy punch	8 mm	25

**Technical specifications:**

- Class IIa
- Sterilization by gamma irradiation
- Single packaging



**Recommended sizes for vulvar biopsy between 2 and 4 mm**

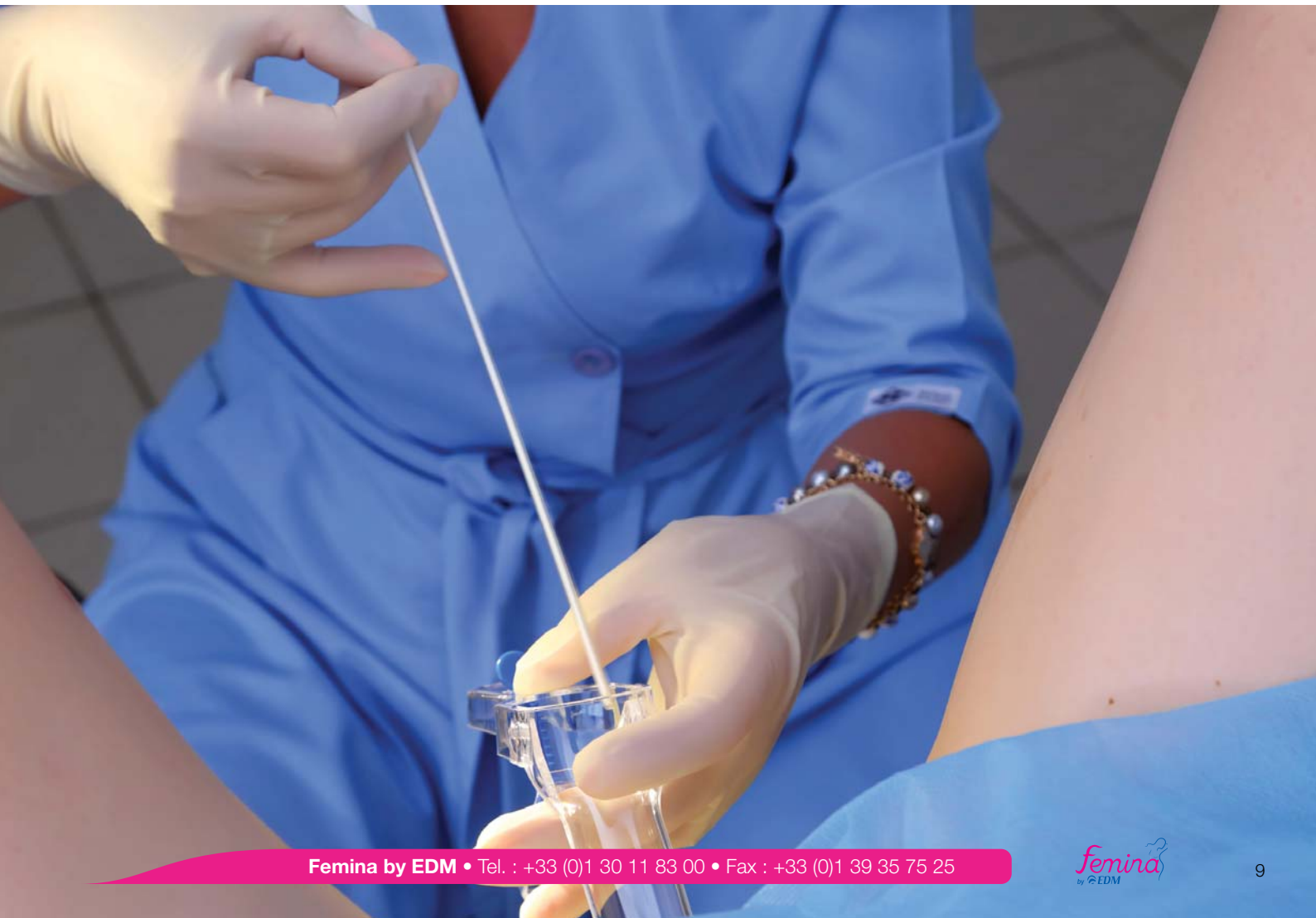



# REUSABLE INSTRUMENTATION

- **Functional exploration**

Cusco and Collin speculum,  
Kogan endocervical speculum ...

- **Stainless steel instrumentation**

Long-handled forceps, Metzenbaum scissors ...

- **Contraception**

Halstead-type forceps, Circeron forceps ...

- **Biopsy and Sampling**

Tischler/Baby Tischler forceps

- **Instrument disinfection**

Aniosyme X3, Disinfection tray



## CUSCO AND COLLIN SPECULUM

- A complete range that adapts to the morphology of all your patients.
- Autoclavable
- Robust



### Foldable Cusco

SKU	Length	Diameter
1404.1010	75 mm	16 mm
1404.1011	80 mm	25 mm
1404.1012	90 mm	30 mm
1404.1013	90 mm	35 mm
1404.1014	90 mm	38 mm
1404.1015	95 mm	40 mm

### Collin

SKU	Length	Diameter
1404.1001	85 mm	25 mm
1404.1002	105 mm	30 mm
1404.1003	110 mm	32 mm
1404.1004	110 mm	35 mm
1404.1005	110 mm	38 mm
1404.1006	110 mm	40 mm
1404.1007	120 mm	30 mm
1404.1008	130 mm	32 mm
1404.1009	130 mm	38 mm

**Technical specifications:**


- Class I non-sterile
- Stainless steel
- Single packaging



## KOGAN ENDOCERVICAL SPECULUM

For visualizing the squamocolumnar junction of the endocervix

- With a ratchet mechanism for secure and effective maintenance of the opening level
- Fine blades
- Reduced patient discomfort

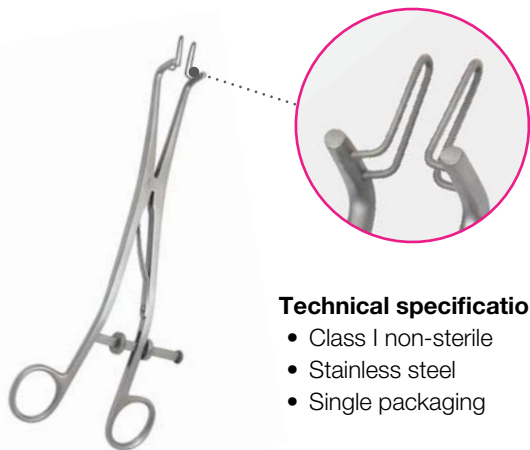

**Technical specifications:**

- Class I non-sterile
- Stainless steel
- Single packaging

SKU	Description	Length	Cdt
1404.0004	Kogan Endocervical Speculum	240 mm	Unité



## DOUBLE HEGAR DILATOR

**For extemporaneous dilation or measurement of the diameter of the cervix**

- Tactile evaluation of the cervix
- 8 sizes available

**INDICATIONS**

- Cervical dilation: pre-abortion, prior to hysteroscopy
- Diagnosis of cervical canal patency (ability to pass a Hegar dilator #8)
- Vaginismus



SKU	Description	Length	Qty
1402.0007	Double Hegar Dilator 5 and 6 mm	200 mm	1
1402.0008	Double Hegar Dilator 11 and 12 mm	205 mm	1
1402.0009	Double Hegar Dilator 13 and 14 mm	205 mm	1
1402.0010	Double Hegar Dilator 19 and 20 mm	205 mm	1

**Technical specifications:**

- Class I non-sterile.
- Stainless steel.

## SILICONE VAGINAL DILATOR SET OF 5

**Designed for independent use by your patients as a complement to physiotherapy or osteopathy sessions**

- Ergonomic and anatomical shape: for easy handling thanks to the external notch
- Made of silicone, soft to the touch
- Slightly flexible dilator for gentle insertion
- Progressive use: 5 different sizes

**INDICATIONS**

- For patients suffering from vaginismus, dyspareunia, vulvodynia, atrophic vaginitis, or Rokitansky syndrome
- Post-vaginoplasty or episiotomy
- During radiotherapy or chemotherapy
- After a difficult menopause



SKU	Couleur	Insertion Length	Total Length	Diameter
1404.0193		90 mm	122 mm	12/14 mm
		100 mm	135 mm	16/18 mm
		110 mm	150 mm	20/22 mm
		120 mm	163 mm	24/26 mm
		130 mm	177 mm	28/30 mm


**Technical specifications:**

- Class I non-sterile
- Medical-grade silicone, latex-free
- Single packaging

Stainless steel instrumentation



Exclusivity

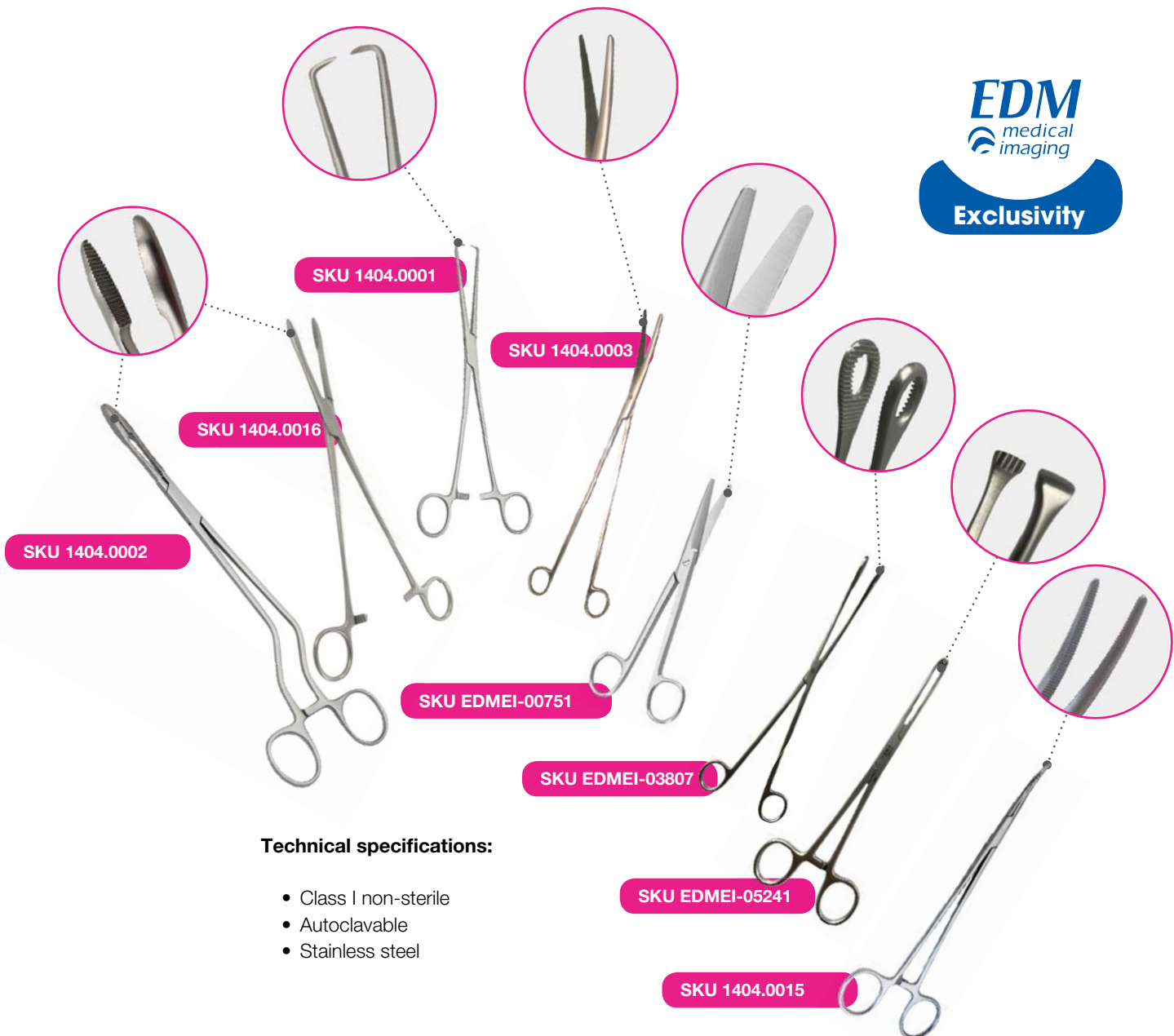

Technical specifications:

- Class I non-sterile
- Autoclavable
- Stainless steel

**LONG-HANDLE FORCEPS (OR CHERON)**

*Straight or curved scissors forceps for grasping dressings and compresses*

- Long handle (250 mm), 2 large oval and serrated jaws
- Ratchet locking system
- Curved version allowing maximum visualization of the cervix

**INDICATIONS**

- Used as a dressing forceps for cervical disinfection or for IUD removal
- Grasping forceps due to its length and flat ends

SKU	Description	Length	Qty
1404.0016	Straight Cheron forceps	250 mm	1 unit
1404.0002	Curved Cheron forceps	250 mm	1 unit

**POZZI FORCEPS**

*Uterine traction to facilitate access to the uterine cavity*

- Facilitates access to the uterine cavity and prevents from damage to the uterus during the introduction of medical instruments

**INDICATIONS**

- Grasping of the cervix and performing a 'traction-fixation' to adjust the uterine flexion angle to 180°.

SKU	Description	Length	Qty
1404.0001	Pozzi forceps	250 mm	1 unit



## METZENBAUM SCISSORS

### Your precision instrument

- Precision instrument for cutting delicate tissue
- Curved and sharp blades

SKU	Description	Length	Qty
1404.0003	Metzenbaum scissors	180 mm	1 unit



INDICATIONS

- Tissue resection: to remove small tumors, polyps, or other benign growths present in the uterus, vagina, or cervix
- Ideal for cutting dressings, bandages
- Recommended for cutting the IUD thread after insertion

## MAYO SCISSORS

### A classic pair of scissors for any gynecology clinic

- Sharp and pointed blade allowing for cutting suture thread close to the wound
- Clean and precise cutting

SKU	Description	Length	Qty
EDMEI-00755	Mayo scissors - curved	170 mm	1 unit
EDMEI-00751	Mayo scissors - Straight	170 mm	1 unit



INDICATIONS

- Cutting and severing suture threads
- Hysterectomy
- Repair of tears or episiotomies
- Excision of lesions or polyps
- During cesarean section

## KELLY FORCEPS

- Slightly curved arms and serrated jaws
- Fenestrated jaws for atraumatic use
- Blunt tip that prevents injuries
- Significant length allowing deep access to the uterus

SKU	Description	Length	Outer diameter of the jaw	Inner diameter of the jaw	Qty
EDMEI-03807	Kelly forceps	320 mm	20 mm	17 mm	1 unit



INDICATIONS

- For late-term abortion

## ALLIS FORCEPS

### Firm tissue grip

- Atraumatic forceps
- Jaws 4 x 5 mm
- Ratchet locking system

SKU	Description	Length	Qty
EDMEI-05241	Allis forceps	190 mm	1 unit



INDICATIONS

- Ovarian surgery
- Hysterectomy
- Repair of tears or episiotomies
- Cervical biopsy

## BENGOLEA FORCEPS

### Hemostatic forceps

- Convenient and easy to use with its long handle (260 mm) and small diameter
- Jaws with claws for effective gripping and holding
- 2 types of models: straight / curved

SKU	Description	Length	Qty
1404.0014	Straight Bengolea forceps	260 mm	1 unit
1404.0015	Curved Bengolea forceps	260 mm	1 unit




INDICATIONS

- For clamping vessels
- Apply traction on the metal suture wire during routine care

**Technical specifications:**

- Class I non-sterile
- Autoclavable
- Stainless steel

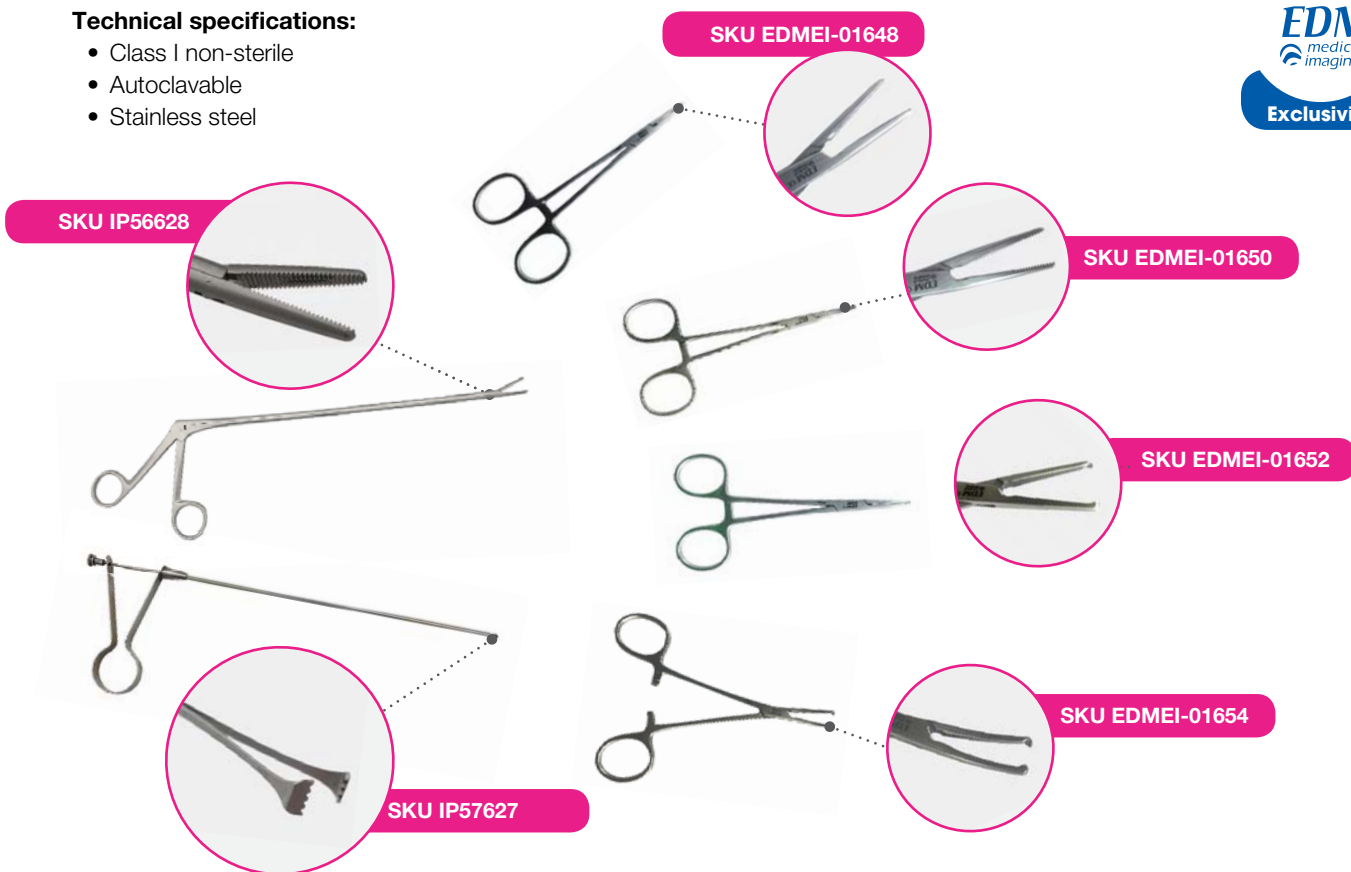

## HALSTEAD-TYPE FORCEPS (HALSTEAD MOSQUITO)

**For the removal of subcutaneous contraceptive implants**

SKU	Description	Length	Qty
EDMEI-01648	Straight forceps without teeth	120 mm	Unité
EDMEI-01650	Curved forceps without teeth	120 mm	Unité
EDMEI-01652	Straight forceps with 1x2 teeth	120 mm	Unité
EDMEI-01654	Curved forceps with 1x2 teeth	120 mm	Unité



- Intended for hemostasis
- Facilitates grasping of devices (dressings, swabs,)
- Removal of subcutaneous contraceptive implant (straight with claws)

## CIRCERON FORCEPS

**For easier removal of the intrauterine device (IUD), especially when the threads are no longer visible or have retracted into the cervix**

- Facilitated introduction into the cervix due to the slimness of this forceps
- Ergonomic: curved handle and clawed ends

SKU	Description	Length	Qty
IP56628	Circeron forceps	280 mm	1 unit

## TERRHUN FORCEPS

**For grasping the intrauterine device (IUD) for insertion or removal**

- Allows for hooking the IUD when the threads are absent
- Fully detachable

SKU	Description	Length	Qty
IP57627	Terrhun forceps	270 mm	1 unit

# Cervical cancer, did you know?



**99%**

of cervical cancers are due to a genital infection caused by oncogenic strains of the human papillomavirus (HPV).



**32,000**

precancerous or cancerous lesions detected every year.<sup>2</sup>



## Cancer screening campaign

Implemented in 2018

17.8 million women aged 25 to 65 affected<sup>3</sup>

90% of cervical cancers could be prevented through screening<sup>3</sup>

## TISCHLER/BABY TISCHLER FORCEPS

- Maximal sampling due to a single tooth on both upper and lower jaws
- Specially designed jaws that minimize the risk of trauma
- Sample fixation with a locking system
- Ergonomic handle



SKU 1404.0005



SKU 1404.0006




- Sampling of tissue from the endocervix, upper part of the cervical canal

SKU	Description	Length	Jaw	Qty
1404.0006	Tischler biopsy forceps	200 mm	3 X 5 mm	1 unit
1404.0005	Baby Tischler biopsy forceps	200 mm	2 X 5 mm	1 unit

## ROTARY TISCHLER/BABY TISCHLER FORCEPS

**Biopsy made easier with these rotary forceps!**

- Sampling of tissue from any angle (up, down, left, right)
- Specially designed jaws that minimize the risk of trauma
- Sample fixation with a locking system
- Ergonomic handle

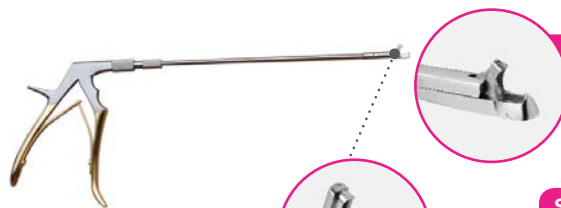

SKU 1404.0010

REF1404.0009



- Cervical biopsy, vaginal biopsy, or even vulvar biopsies
- Excision of small lesions

SKU	Description	Length	Jaw	Qty
1404.0009	Rotary Tischler biopsy forceps	200 mm	3 X 7 mm	1 unit
1404.0010	Rotary Baby Tischler biopsy forceps	200 mm	2 X 5 mm	1 unit


1. <https://www.euro.who.int/fr/health-topics/noncommunicable-diseases/cancer/news/news/2012/2/early-detection-of-common-cancers/cervical-cancer>  
 2. Dépliant - Dépistage du cancer du col de l'utérus. Guide pratique DEPCCU21. [www.e-cancer.fr](http://www.e-cancer.fr)  
 3. <https://www.e-cancer.fr/Comprendre-prevenir-depister/Se-faire-depister/Depistage-du-cancer-du-col-de-l-uterus>



## KEVORKIAN FORCEPS

- Narrow jaws facilitating insertion and sampling in the cervical canal
- Sample fixation with a locking system

SKU	Description	Length	Jaw	Qty
1404.0007	Kevorkian biopsy forceps	245 mm	3 x 10 mm	1 unit

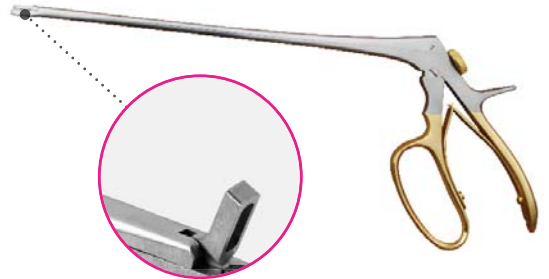

## TISCHLER-KEVORKIAN FORCEPS

- Cervical biopsy punch
- Sample fixation with a locking system

SKU	Description	Length	Jaw	Qty
1404.0008	Tischler-Kevorkian forceps	230 mm	3 x 8 mm	1 unit



## NOVAK TYPE BIOPSY CANNULA

- Very fine
- Equipped with small teeth to facilitate grasping of IUD threads



- Endometrial biopsy
- Aid tool for the removal of the intrauterine device in case of migration into the uterine cavity

SKU	Description	Length	Diameter	Qty
1404.0011	Novak cannula with incorporated reservoir	230 mm	3 mm	1 unit



## KEVORKIAN CURETTE OR UTERINE CURETTE


**Used for uterine biopsies and endometrial curettage by scraping or aspiration**

- Edges are curved inwards and slightly blunt
- Ratchet design to maintain consistent opening level during the procedure
- Ultra-thin valves to minimize patient discomfort
- Without basket



- Used in colposcopy to perform samplings
- Removal of uterine polyps
- Treatment of excessive bleeding
- Uterine cleaning after an incomplete miscarriage

SKU	Description	Length	Diameter	Qty
1404.0012	Kevorkian curette or uterine curette	300 mm	4 mm	1 unit

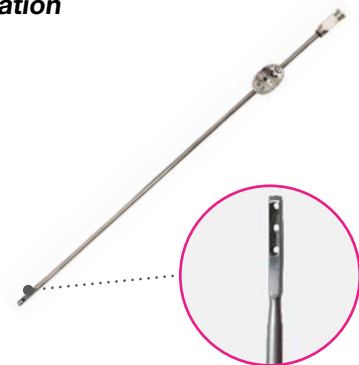

All our range of reusable instruments must be cleaned, disinfected, and sterilized before use.

## 1 Pre-disinfection (immersion)

This step involves removing any organic residues, such as blood or tissues. It includes soaking the instruments in an enzymatic solution to break down proteins and organic matter. Be sure to observe the soaking time.

### Our selection

ALDEHYDE-FREE



## ANIOSYME X3

**Soaking time: 5 min for X3 and 15 min for DD1**

- Cleaning and pre-disinfection of medical and surgical instruments, medical devices, and endoscopic equipment
- Used at a dilution of 0.5% - Neutral pH

### Microbiological properties

**BACTERICIDE:** EN 1040, EN 13727, EN 14561, MDR according to EN 13727: active in 5 minutes / 0,5%.  
**LEVURICIDE:** EN 1275, EN 13624, EN 14562: active in 5 minutes / 0,5%.  
**VIRUCIDE:** EN 14476: active against HIV-1, PRV (model virus HBV), Herpesvirus, BVDV (virus modèle HCV), Vaccinia virus: active in 5 minutes / 0.5%.

SKU	Description	Qty
AN2633095	1 liter with built-in dispenser	1 unit
AN2633036	5 liters with built-in dispenser	1 unit

### Add-On Product



## DISINFECTATION TRAY INSTRUBAC

- Soaking tray with a useful volume of 5 liters or 2 liters
- Includes a solid lid without slots and a polypropylene basket for collecting and transporting soiled instruments.

SKU	Description	Dimensions	Qty
404002NM	Soaking tray 10 L Effective volume 5 L	Overall dimensions with lid: 43 cm (L) x 29 cm (W) x 12.5 cm (H) Dimensions of the tray: 36.5 cm (L) x 25.5 cm (W) x 9 cm (H)	1 unit
AN404038	Soaking tray 5 L Effective volume 2 L	Overall dimensions with lid: 29.5 cm (L) x 20 cm (W) x 13 cm (H) Dimensions of the tray: 23 cm (L) x 17 cm (W) x 9 cm (H)	1 unit

## 2 Rinse

After pre-disinfection, the instruments should be thoroughly rinsed with water to remove any residue or particles detached during the cleaning process.



Disinfection of instruments

3

Drying

The instruments should be thoroughly air-dried on a clean work surface to prevent any subsequent contamination, before autoclave sterilization or cold disinfection.

4

Disinfection or Sterilization

Actual disinfection involves using an appropriate disinfectant to kill any remaining microorganisms on the instruments. It can be done cold with a suitable disinfectant or by steam autoclave sterilization (type class B) after packaging in specific pouches or boxes.

In gynecology, commonly used disinfectants are chemical products like glutaraldehyde, hydrogen peroxide, ethanol, and isopropanol.

Our selection

ANIOXYDE 1000


Immersion time: 5 min for intermediate and high-level disinfection

- High-level disinfection of medical devices
- Extemporaneous production of peracetic acid (without acetic acid)
- Glutaraldehyde-free
- Effectiveness of the bath verified by test strip

Microbiological properties

**BACTÉRICIDE:** EN 1040, EN 14561, EN13727 active in 5 min.  
**LEVURICIDE / FONGICIDE:** EN 14562, EN 13624, EN 1275 active in 5 min.  
**VIRUCIDE:** EN 14476 active in 5 minutes.  
**MYCOBACTÉRICIDE:** EN 14348, EN14563: active in 5 min.  
**SPORICIDE:** EN 14347, EN13704: active in 5 min.  
*Retrouvez l'intégralité des propriétés microbiologiques sur la fiche technique disponible sur simple demande.*

ALDEHYDE-FREE

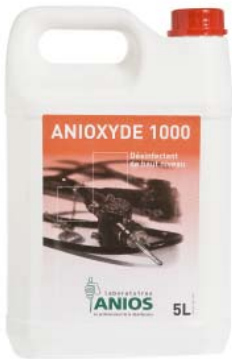

SKU	Description	Qty
AN1081299	5L canister of Anioxyde 1000	1 unit
AN100194	Control strips for Anioxyde 1000	50

5

Storage

Disinfected instruments should be thoroughly dried and stored in a clean, dry place until their next use.

## SELF-SEALING AUTOCLAVE STERILIZATION POUCH

### 2 in 1, pouch with sterilization indicator

- Perfect sealing of the pouch thanks to the two layers of adhesive on the flap
- Integrated sterilization indicator



- Sterilization for moist heat or ethylene oxide.

#### Technical specifications:

- Compliant with European and international standards: MDR (EU) 2017/745 - EN ISO 11607-1 / 11607-2 - EN 868-5

SKU	Description	Qty
EDME5011	Sterilization pouch 90 x 250 mm	9 boxes – box of 200 sachets
EDME5012	Sterilization pouch 140 x 260 mm	6 boxes – box of 200 sachets
EDME5015	Sterilization pouch 190 x 330 mm	6 boxes – box of 200 sachets

## STANDARD STERILIZATION INDICATOR

### Check the effectiveness of steam sterilization

- Ensure optimal sterilization results with autoclave and reliable traceability
- Combination of sterilization parameters: temperature, duration, saturated steam
- Exposure to a sterilization cycle changes the color of the indicator from yellow to black after
- 15 minutes at 121°C or 5 minutes and 30 seconds at 134°C



#### Technical specifications:

- Compliant with standards EN867-5 and ISO11140-1 and 4

SKU	Description	Qty
EDME5110	Sterilization indicator	250

## PRION TEST

### Sterilization cycle control test

- Ensures sterilization under optimal conditions for complete prion destruction
- Clear evidence of exposure to a sterilization cycle by changing color from yellow to black under the following conditions: 18 minutes at 134°C



SKU	Description	Qty
EDME5116	Prion sterilization indicator	250

## HELIX TEST

### For hollow item monitoring and ensuring optimal sterilization

- Ensures proper functioning of steam sterilizers by verifying steam penetration into hollow items
- Reacts after 3.5 minutes at a temperature of 134°C or 15 minutes at 121°C

#### Technical specifications:

- Can only be used in autoclaves with a class B cycle (fractionated vacuum)
- Compliant with standards EN 867-5, EN ISO 11140-1 (class 6) & 11140-4



SKU	Description	Qty
EDME5115	Helix test	100


# WOMEN'S IMAGING

- **Hysterosalpingography**  
Catheter, probe, cannula
- **Imaging Consumables**  
Ultrasound gel, probe cover
- **Probe Disinfection**  
Chemical and automatic disinfection
- **Radiology and X-ray**  
Cartridge, Thermal paper, printer

## BALLOON CATHETER

### Sterile with balloon


- Easy and simplified handling
- Sealing ensured by the silicone balloon's performance
- Atraumatic insertion: the insertion sheath eliminates the need for a tenaculum.

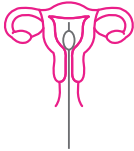

- Designed for use in sonohysterography and hysterosalpingography, the catheter can be used with saline-based media according to the procedure's requirements.

#### SKU 61-5005

The kit includes an insertion sheath, a syringe for inflating the balloon, and a clamp for blocking the injection pathway. The protective sheath is made of polyvinyl chloride, and the insertion sheath is made of polypropylene.

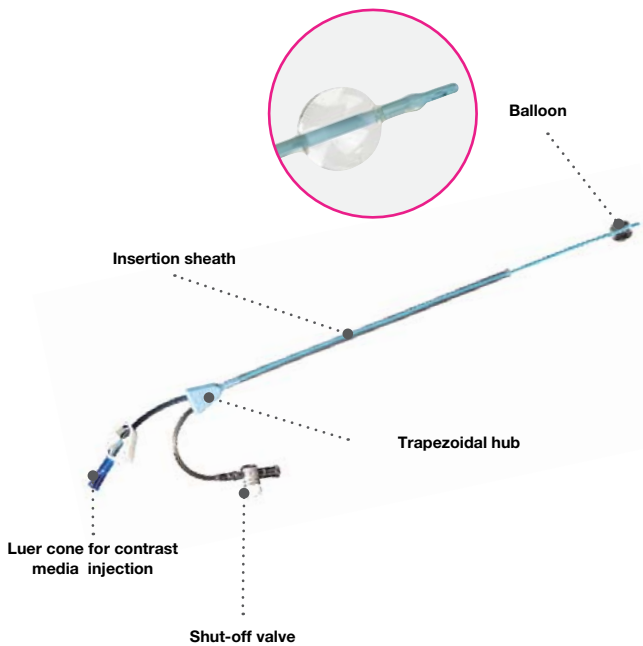

SKU	Catheter diameter	Syringe volume	Length	Qty
61-5005	5Fr	1,5 cc	300 mm	10
61-5007	7Fr	3 cc	300 mm	10

#### Technical specifications:

- Class I sterile
- Ethylene oxide sterilization
- Single use
- Polyurethane catheter
- Latex-free
- Available in 5 and 7 Fr

## DOUAY HYSTEROSCOPE COMPLETE SET REUSABLE

- This kit includes:  
A stainless steel box measuring 30 x 15 x 7 cm, a Douay uterine cannula with 3 cones of sizes 20 mm, 25 mm, and 35 mm, a 20 cc Douay syringe, a 24 cm stainless steel Pozzi clamp, a 24 cm stainless steel dressing forceps, a Collin speculum 30 x 90 mm.



#### Technical specifications:

- Class I
- Reusable

SKU	Description	Qty
HYS1010	Reusable Complete Douay Hysteroscope	1 unit

## SONY COLOR PAPER KITS



SKU	Kit contents	Printer compatibility	Qty
<b>UPC-21S</b>	Color Printing Kit A7: 240 sheets of 100 x 90 mm, 3 ink ribbons	UP-20, UP-25MD, UP-D25MD, UP-21MD, UPD-23MD	1 unit
<b>UPC-21L</b>	Color Printing Kit A6: 200 sheets of 140 x 100 mm, 4 ink ribbons		1 unit
<b>UPCR80MD</b>	Laminated Color Printing Kit A4: 100 sheets (50 sheets x 2 rolls), 2 ink ribbons	UPDR80MD	1 unit

## SONY UP-X898MD

**High-quality Black & White Printing - A6**



### A6

Noir & blanc

### Consumables

UPP-110HG

UPP-110HD

UPP-110S



INDICATIONS

- Ultrasound
- Mobile gantry
- Echo-endoscopy

The Sony UP-X898MD thermal printer is the ideal solution for digital ultrasound systems.

- High-quality printing with high resolution, approximately 325 dpi, and grayscale levels of up to 256 steps
- Fast printing speed of approximately 2 seconds
- Selectable image format: 4:3 or 1:1
- Compact design with reduced footprint



## THERMAL PAPERS FOR PRINTING ULTRASOUND IMAGES

- Compatible with Mitsubishi printers
- A6 format
- Black and white printing



### MITSUBISHI PAPERS

- Compatible with Mitsubishi printers
- Black and white printing

110 mm (l)  
x 18 m (L)




SKU	Description	Format du papier	Compatibility	Qty
<b>K61B</b>	Matte paper	110 mm (W) x 20 m	P93E	Per 4 rolls
<b>K65HM</b>	High-density paper	110 mm (W) x 20 m	P93E	Per 4 rolls
<b>K91HG</b>	High-gloss paper	110 mm (W) x 20 m	P93E	Per 4 rolls
<b>K95HG</b>	High-gloss paper	110 mm (W) x 18 m	P95DW	Per 5 rolls

### SONY PAPERS

- Compatible with Sony printers
- A4, A6, A7 format
- Black and white printing

110 mm (l)  
x 18 m (L)




SKU	Description	Compatibility	Qty
<b>UPP110 HG</b>	Glossy paper	UP-897MD, UP-D897, UP-X898MD	10
<b>UPP110 HD</b>	High-density paper	UP-897MD, UP-D897, UP-X898MD	10
<b>UPP110 S</b>	Standard paper	UP-897MD, UP-D897, UP-X898MD	10

210 mm (l)  
x 25 m (L)



SKU	Description	Compatibility	Qty
<b>UPP210 HD</b>	High-density paper	P-910AD, UP-930AD, UP-960AD, UP-970AD, UP-980AD, UP990AD	5
<b>UPP210 SE</b>	Standard paper	UP-991AD, 990AD, 980AD, 971AD, 970AD, 960AD, 930AD, 910AD	5

84 mm (l)  
x 12,5 m (L)



SKU	Description	Compatibility	Qty
<b>UPP84 HG</b>	Glossy paper	UP-D711MD	10
<b>UPP84 S</b>	Standard paper	UP-D711MD	10



# How to correctly choose your gel and probe cover



## NON-CRITICAL

## SEMI-CRITICAL

## CRITICAL

### PROCEDURE

Surface ultrasound on intact skin

Endocavitary ultrasound  
Surface ultrasound on wounds and burns

Intraoperative ultrasound  
Anesthesia, CVC

Probe in contact with intact skin

Probe in contact with mucosa or injured skin

Probe in contact with sterile tissue or cavities or vascular system

### EXAMPLES IN ULTRASOUND



### GEL

Non-sterile gel in a 250 ml bottle or a 20 ml sachet

Sterile gel in a 20 ml sachet

Sterile gel in a 20 ml sachet

### PROBE PROTECTION

Not required

Sterile protection

Sterile protection

# Did you know?



The use of sterile gel is highly recommended for all semi-critical and critical US procedures, such as transducer contact with mucous membranes, i.e., all endo-cavity US, contact with any body fluids, i.e., all major and minor US guided interventional procedures, and when scanning infected or broken skin and wounds.

Nyhsen CM, Humphreys H, Koerner RJ, Grenier N, Brady A, Sidhu P, Nicolau C, Mostbeck G, D'Onofrio M, Gangi A, Claudon M. Infection prevention and control in ultrasound - best practice recommendations from the European Society of Radiology Ultrasound Working Group. *Insights Imaging*. 2017 Dec;8(6):523-535. doi: 10.1007/s13244-017-0580-3. Epub 2017 Nov 27. PMID: 29181694; PMCID: PMC5707224.

Specializing in ultrasound gels for over 30 years, EDM Imaging has developed its own range of non-sterile gels to meet the various needs of practitioners. Compliant with the new Regulation 2017/745.

## NON-STERILE ULTRASOUND GEL

The SONECHO, PREMIUM, POLARIS II ranges from the EDM group are ultrasound gels manufactured entirely in France. We are committed to providing you with the highest possible product quality, for a healthy gel without artifacts.



## NON-STERILE ULTRASOUND GEL

- High viscosity for excellent ultrasound transmission
- Non-irritating, pH neutral
- Fragrance-free
- Preserves probe integrity
- Our gels are easy to apply and clean, without staining

SKU	Range	Color	Specification	Volume	Qty
PC250	Premium Gel	Clear	-	250 ml bottle	24 or 48
CG5000	Premium Gel	Clear	-	5 L flexible cubitainer with pouring spout	1 unit
PC250/24B	Sonecho gel	Blue	-	250 ml bottle	24
CG5000B	Sonecho gel	Blue	-	5 L flexible cubitainer with pouring spout	1 unit
PC250/24SP	Polaris gel	Clear	Paraben-free, hypoallergenic	250 ml bottle	24
CG5000SP	Polaris gel	Clear	Paraben-free, hypoallergenic	5 L flexible cubitainer with pouring spout	1 unit
KG250/40	Krystal gel	Clear	Paraben-free, hypoallergenic	250 ml bottle	40
PA60	Parker Aquasonic gel	Blue	-	60 gr Tube	12
PA250	Parker Aquasonic gel	Blue	-	250 ml bottle	12
PA1000	Parker Aquasonic gel	Blue	-	1 L bottle	1 unit
PA5000	Parker Aquasonic gel	Blue	-	5 L cubitainer	1 unit
PA251	Parker Aquasonic gel	Clear	-	250 ml bottle	12
PA5010	Parker Aquasonic gel	Clear	-	5 L cubitainer	1 unit

## STERILE ULTRASOUND GEL

- Double sterile packaging
- Narrow and precise distribution
- Completely aqueous, no risk of staining

SKU	Description	Capacity	Qty
PC020	Krystal by EDM	20 ml	48
PA020	Parker Aquasonic 100	20 ml	48

Mandatory use as per Circular No. 960479 of the Ministry of Labor and Social Affairs dated 06/02/1996.



### Technical specifications:

- Class I
- Sterilization by Gamma rays
- 20 ml sachet

## Hygiene and disinfection recommendations for endocavitary ultrasound procedures

(protocol based on SF2H recommendations, March 2019)

### STEP 1

- Hand hygiene / Hydroalcoholic friction
- BEFORE and AFTER the examination
- Application of sterile gel on the probe (20 ml per sachet)
- Placing of the sheath (CE)



### STEP 2

Wearing single-use GLOVES



### STEP 3

At the end of the procedure:

- Gently remove the used cover
- Dispose of the gloves
- Perform hand hygiene, put on non-sterile single-use gloves, and wipe off the remaining gel with a non-abrasive absorbent material



### STEP 4

Clean the probe by wiping it with a detergent wipe.

- Proceed with an intermediate-level disinfection.





Probe cover



Non-sterile probe cover for endocavitary

EDM ECHO COVER™

Latex

Endocavitary ultrasound probe covers. The most economical solution for your consultations.



SKU	Dimensions	Description	Folding	Sterile	Qty
600249	3,5 x 20 cm	Non-lubricated EDM	Rolled	NO	144
599249	3,5 x 20 cm	Lubricated EDM Protection	Rolled	NO	144
650249	3 x 30 cm	Non-lubricated EDM Protection	Rolled	NO	100

Polyisoprene - latex-free

SKU	Dimensions	Description	Folding	Sterile	Qty
670249	3 x 30 cm	Non-lubricated EDM Protection	Rolled	NO	100



Sterile superficial probe cover

KRYSTAL

Polyethylene

Flexible polyethylene protections, offering an excellent value for money for sterile procedures, especially the more invasive ones.



SKU	Dimensions	Description	Folding	Sterile	Qty
YPE1830	18 x 30 cm		Telescopic	Sterile EO*1	100
YPE1880	18 x 80 cm	3D end covers, supplied with 2 fixing elastics and 1 drape	Telescopic	Sterile EO*1	100
YPE18125	18 x 125 cm		Telescopic	Sterile EO*1	50
YPE1530WG	15 x 30 cm		Telescopic	Sterile EO*1	25
YPE1580WG	15 x 80 cm	Echo Cover™ Protection Kit (1 cover, 1 drape, 1 sterile gel sachet, 2 fixing elastics)	Telescopic	Sterile EO*1	25
YPE15125WG	15 x 125 cm		Telescopic	Sterile EO*1	25
YPE15240WG	15 x 240 cm		Telescopic	Sterile EO*1	25

KRYSTAL COVER

Kit - Polyurethane

- Excellent value for money - Elasticity and strength - The kit includes: 1 sterile cover, 1 sterile drape, 2 sterile elastics, and 1 individual sterile gel sachet 20 g

STERILE



SKU	Dimensions	Description	Folding	Sterile	Qty
PU1590WG	15 x 90 cm	Protection kit	Telescopic	Yes	25
PU1560WG	15 x 60 cm	Protection kit	Telescopic	Yes	25
PU1530WG	15 x 30 cm	Protection kit	Telescopic	Yes	25
PU15240WG	15 x 240 cm	Protection kit	Telescopic	Yes	25
PU15150WG	15 x 150 cm	Protection kit	Telescopic	Yes	25
PU15125WG	15 x 125 cm	Protection kit	Telescopic	Yes	25

SKU	Dimensions	Description	Folding	Sterile	Qty
PUI1560WG	15 x 80 cm	Protection Kit with Insert	Telescopic	Yes	25
PUI15125WG	15 x 125 cm	Protection Kit with Insert	Telescopic	Yes	25

# Did you know ?




**SF2H Recommendations, March 2019**

Prevention of infectious risk associated with endocavitary ultrasound procedures

« Intermediate-level disinfection in accordance with manufacturers' recommendations using one of the following methods:

- Automated intermediate or higher-level processes
- Disinfectant wipes
- Any other validated intermediate-level disinfection process»

## Levels of disinfection required according to the Spaulding Classification

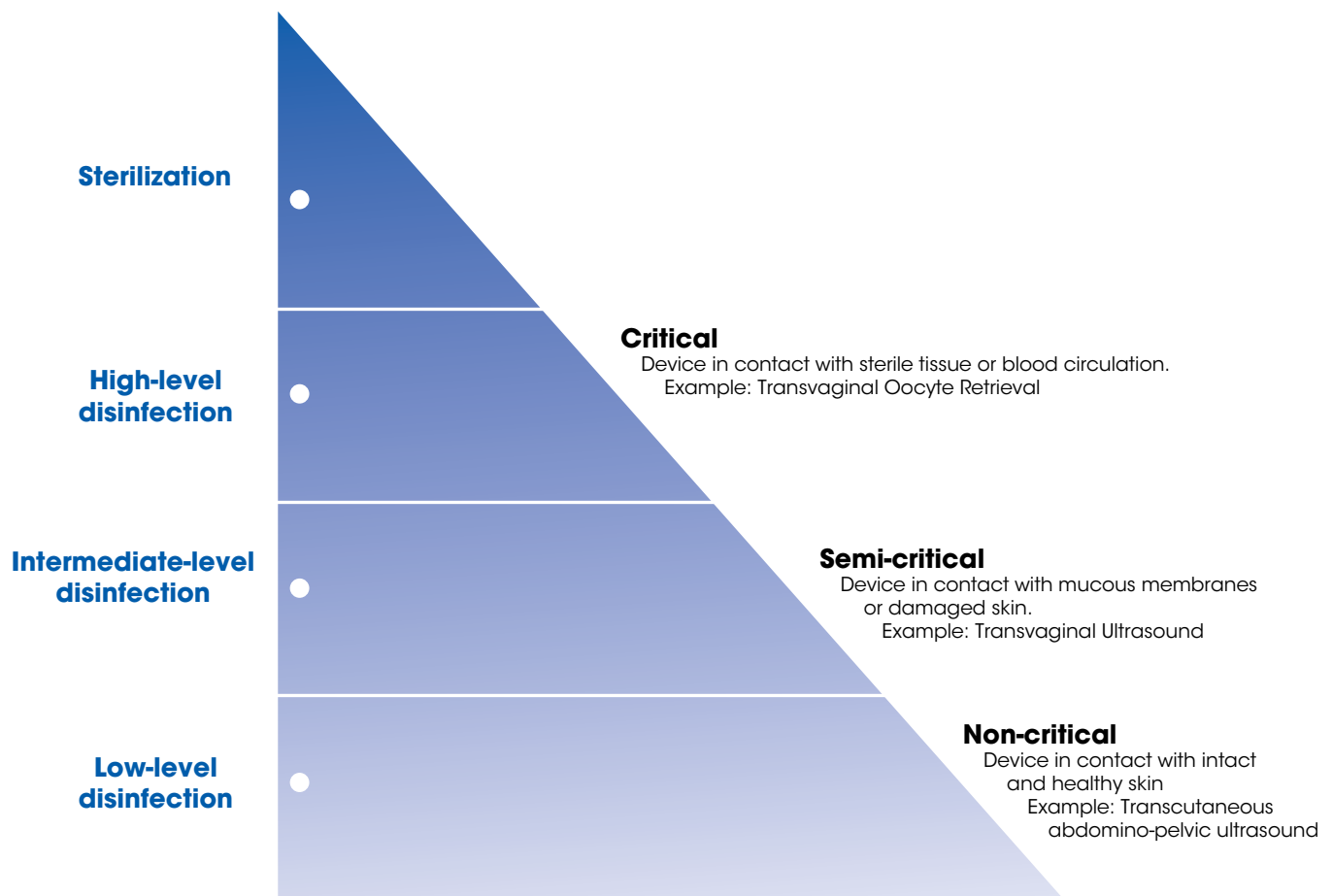

Low-level disinfection



ALCOHOL-FREE

**SEPTIWIPES® / SEPTISPRAY®**

Alcohol-free spray & impregnated wipe

Full activity in 2 minutes for wipes and 5 minutes for spray

- For cleaning and disinfecting ultrasound probes
- Aldehyde-free, hypoallergenic

Microbiological properties:

	Spray	Impregnated wipe
<b>BACTERICIDE:</b>	EN 13727, EN 1276, EN 13697, EN 14561, NF T 72-281: active in 5 min.	EN 16615, EN 1275, EN 13624 active in 1 min.
<b>YEASTS AND MOLDS:</b>	EN 13624, EN 1650, EN 13697, EN 14562, NF T 72-281 : active in 5 min.	EN 1275, EN 13624
<b>VIRUCIDE:</b>	EN 14476 active in 5 min.	EN 14476 active in 2 min.
<b>SPORICIDE:</b>	NF T 72-281 active in 5 min.	

Find all the microbiological properties on the technical sheet available on request

Septiwipes® wipes are certified for probe disinfection ultrasound by:

SKU	Description	Qty
SP2000	Dispenser box of 100 wipes, 200 x 140 mm	1 unit
SP750	750 ml spray bottle	1 unit

Intermediate-level disinfection



**ECHO CLEAN®**

Total activity in 5 minutes

- For the cleaning and disinfection of ultrasound probes
- Aldehyde-free, hypoallergenic
- Active against Papillomavirus (human strains HPV-18 and HPV-16)

Microbiological properties:

<b>BACTERICIDE:</b>	- EN 13727: active in 1 min. - EN 16615 (with mechanical action under dirty conditions): Pseudomonas aeruginosa, Staphylococcus aureus, Enterococcus hirae: active in 5 min.
<b>YEASTS/MOLDS:</b>	- EN 13624: active in 1 min. - EN 1650, EN 13697, EN 14562, EN 16615: active in 5 min.
<b>VIRUCIDE:</b>	- EN 14476: active in 5 min. Protocole Ph. D Meyers' HPV-18 (>4log), HPV-16 (3,5log) <sup>2</sup> 5 min
<b>SPORICIDE:</b>	- EN 13704 active in 15 min.

Find the full microbiological properties on the technical data sheet available upon request.

1 Protocole adapté de la norme EN 16615 par le PennState College of Medicine Department of Microbiology and Immunology Pennsylvania, Ph. D Craig Meyers. 2 Protocole spécifique aux lingettes

Echo Clean is certified for the disinfection of ultrasound probes by:

SKU	Description	Qty
EC2000	Dispenser box of 120 wipes, 130 x 190 mm	1 unit
EC1000	1 L spray bottle	1 unit

High-level disinfection



**VIROWIPES**

Total activity in 15 minutes

- Cleaning and disinfection of medical devices and ultrasound probes
- Without CMR (Carcinogens, Mutagens, toxic to Reproduction), without aldehyde, without quaternary ammonium
- Complete sporicidal activity

Microbiological properties:

<b>BACTERICIDE:</b>	EN 13727, EN 1276, EN 13697: active in 5 min.
<b>YEASTS/MOLDS:</b>	EN 13624, EN 1650, EN 13697: active in 5 min.
<b>VIRUCIDE:</b>	EN 144476: Norovirus Murin: active in 5 minutes. EN 14476: Adénovirus, Rotavirus et Poliovirus: active in 10 minutes.
<b>TUBERCULOCIDE:</b>	EN 14348: active in 10 min.
<b>SPORICIDE:</b>	EN 13704 (dont Clostridium Difficile): active in 15 min. EN 17126: bacillus subtilis : active in 10 min. EN 17126: bacillus cereus : active in 15 min.

Find the full microbiological properties on the technical data sheet available upon request.

SKU	Description	Qty
WP2000	Box of 150 Virowipes cleaning/disinfecting wipes	1 unit

# CONSUMABLE

- **Sheet and cotton**  
Examination sheet
- **Protection and dressing**  
Examination gloves, tunic and gown, cap, mask
- **Care**  
Dressings, cotton, dressing fields, scalpel, needle, syringe, dressing set, adhesive tape, needle collector
- **Maintenance**  
Surface disinfection
- **Hygiene**  
Hand hygiene, intimate hygiene



# Hygiene and disinfection solutions for your medical practice

**NuClean**  
Spray and wipes for disinfecting surfaces of medical devices  
P. 44

**NuClean SH Surfa'safe**  
Spray for disinfecting high surfaces  
P. 44

**Examination drapes**  
P. 34

**Absorbent towels/tissues**  
P. 34

**Viro'wipes**  
Wipes for high-level disinfection of probes and medical devices  
P. 31

**EchoClean**  
Spray and wipes for intermediate-level disinfection of probes and medical devices  
P. 31

**Septispray & Septiwipes**  
Spray and wipes for low-level disinfection of probes and medical devices  
P. 31

**Sensitex**  
Sterile & non-sterile latex gloves  
P. 37

**Aniosgel 800**  
Hand sanitizer  
P. 45

- Disinfection of probes and medical devices
- Hand hygiene
- Environmental hygiene
- Surface hygiene



## EXAM BED COVER

Soft and durable



### SMOOTH IN WADDING

Ensure protection for all types of medical examinations

SKU	Description	Qty
516	White cotton sheets 2x18 g/m <sup>2</sup> , 150 sheets, 50 cm x 35 cm	12
522	Recycled cotton sheets 2x18 g/m <sup>2</sup> , 150 sheets, 50 cm x 35 cm	12



### EMBOSSED

High resistance and additional comfort for the patient

SKU	Description	Qty
523	Embossed white cotton sheets 2x18 g/m <sup>2</sup> , 135 sheets, 50 cm x 35 cm	9



### PLASTICIZED

A resistant and waterproof plastic-coated side to protect your couch

- Made in France
- Durable: 36 g/m<sup>2</sup>
- Convenient: 1 unit roll is pre-cut
- Comfortable

SKU	Description	Qty
519	Blue plastic-coated cotton sheets 36 g/m <sup>2</sup> , 150 sheets, 50 cm x 38 cm	6



## HAND TOWELS & TOILET PAPER

SKU	Description	Qty
708	Interfolded white 2-ply paper towels (packaged in 30 packs of 125) unfolded 22 x 35 cm, folded 22 x 9 cm	Unité
709	Hand towels 25 x 38 cm, side dispensing	24
710	Hand towels 22 x 38 cm, center dispensing	12
S401023	Toilet paper, 200 sheets/roll	96

## SURGICAL GOWN

- Hypoallergenic
- Excellent comfort
- Short sleeves, single closure with drawstring at the waist
- Non-transparent, appreciated by the patient
- One size fits all

**Technical specifications:**

- Class I non-sterile
- Individual packaging
- Non-woven polypropylene



SKU	Dimensions	Qty
PFS-01	Operated profile Shirt, adult, navy blue	10

## VISITOR GOWN

- Neckline: Tied with 2 laces
- Sleeve: Elastic cuffs

**Technical specifications:**

- Class I non-sterile
- Individual packaging
- Non-woven polypropylene

SKU	Description	Qty
161265/100	Visitor white gown 25g/m <sup>2</sup> - 110 X 130 cm	100
161266/100	Visitor blue gown 25g/m <sup>2</sup> - 110 X 130 cm	100



## NON - STERILE SCRUB WITH 2 POCKETS

- This scrub consists of a short-sleeved tunic and trousers
- Hypoallergenic
- Trousers adjustable with a tie
- V-neck
- Available in 5 sizes

**Technical specifications:**

- Class I non-sterile
- Individual packaging
- Non-woven polypropylene
- Weight: 40 g

SKU	Description	Dimensions in cm Pants length / Chest circum	Size	Qty
PFT-03	Block Profile Blue Scrub	96 / 55	S	50
PFT-04	Block Profile Blue Scrub	101 / 61	M	50
PFT-05	Block Profile Blue Scrub	108,5 / 63,5	L	50
PFT-06	Block Profile Blue Scrub	114 / 72	XL	50
PFT-07	Block Profile Blue Scrub	118 / 79	XXL	50



# How to properly choose your gloves



	VINYLE	LATEX	NITRILE
<b>Use</b>	Short-duration procedures involving neither blood nor bodily fluids.	Use for precise actions requiring sensitivity and fine dexterity.	Use for precise actions requiring sensitivity and dexterity, with the advantage of being less allergenic and more resistant to chemicals.
<b>Protective barrier</b>	★★★★☆	★★★★★★	★★★★★★
<b>Sturdiness</b>	★★★★☆	★★★★★★	★★★★★★
<b>Elasticity</b>	★★★★☆	★★★★★★	★★★★★★
<b>Resistance to puncture</b>	★★★★☆	★★★★★★	★★★★★★
<b>Chemical resistance</b>	★★★★☆	★★★★★★	★★★★★★
<b>Fit and comfort</b>	★★★★☆	★★★★★★	★★★★★★
<b>Tactile sensitivity</b>	Inferior	Excellent	Adequate Less flexible than latex
<b>User Risks</b>	Low because there is no latex protein but presence of chemical ingredients. Ideal for individuals sensitive to allergies.	Associated with the presence of latex, it can cause allergic reactions. Risk of allergy types I and IV.	Low because there is no latex protein but presence of chemical ingredients. Risk of allergy type IV.
<b>Indication</b>	They can be used in medical environments with low risk, such as handling non-hazardous materials or administering medications.	Becoming less common due to allergy concerns. They are sometimes replaced by nitrile gloves.	Healthcare professionals in contact with chemicals, patients allergic to latex, and situations requiring high resistance.

## NITRILE GLOVE

### Non-powdered and non-sterile

- Hypoallergenic
- Anti-slip system: textured
- Ambidextrous - blue - rolled cuff

SKU	Description	Dimensions	Qty
360S	Nitrile examination gloves, latex-free, powder-free - Small 6-7	240 x 80 mm	100
360M	Nitrile examination gloves, latex-free, powder-free - Medium 7-8	240 x 95 mm	100
360L	Nitrile examination gloves, latex-free, powder-free - Large 8-9	240 x 105 mm	100



## LATEX GLOVE

### Powdered and non-powdered, sterile and non-sterile

- Non-slip texture
- Ambidextrous
- Powdering: cornstarch
- Color: milky white - matte appearance

SKU non-powdered	SKU Powdered	Steriles	Dimensions	Taille
340S	340S/P	NO	240 X 80 mm	Small (6/7)
340M	340M/P		240 X 95 mm	Medium (7/8)
340L	340L/P		240 X 110 mm	Large (8/9)
-	3524	YES	240 X 80 mm	Small (6/7)
-	3525		240 X 95 mm	Medium (7/8)
-	3526		240 X 110 mm	Large (8/9)



## VINYL GLOVE

### Powdered and non-powdered, non-sterile

- Ambidextrous
- Color: Transparent
- Powdering: cornstarch

SKU non-powdered	SKU Powdered	Dimensions	Size
3501	3501/P	240 x 80 mm	Small (6/7)
3502	3502/P	240 x 95 mm	Medium (7/8)
3503	3503/P	240 x 110 mm	Large (8/9)



## LEGUEU TYPE 1 OR 2 FINGERS

### Made of polyethylene

- Legueu type, opaque

SKU	Description	Qty
DP422	Non sterile - 2 fingers	100
DS426	Non sterile - 2 fingers, individually packaged	100





## MASK

SKU	Description	Qty
<b>MS301</b>	FFP2 Single-use Mask	20
<b>PLM-01</b>	Pink Surgical Mask - Type II - Non-sterile - 3 folds	50
<b>PLM-02</b>	Blue Surgical Mask - Type II - Non-sterile - 3 folds	50

## GAUZE SWAB

- Does not fray thanks to a surgical fold with tucked edges
- Good absorption capacity
- Woven from hydrophilic cotton yarn, made with 17 threads
- Softness during use



### Technical specifications:

- Class I
- Non-sterile
- 17 threads, 23 g/m<sup>2</sup>
- 8 folds
- Standard: NF EN 14079



- Cleaning and disinfection operation of the wound
- Drying of a lesion thanks to its high hydrophilicity
- Constituent element of the dressing

### Non Sterile

SKU	Dimensions	Qty
<b>18500</b>	5 x 5 cm	100
<b>18503</b>	7,5 x 7,5 cm	100
<b>18506</b>	10 x 10 cm	100

### Sterile

SKU	Dimensions	Qty
<b>13600</b>	5 x 5 cm	Box of 100 (50 pouches of 2 swabs)
<b>13601</b>	7,5 x 7,5 cm	Box of 80 (40 pouches of 2 swabs)
<b>13602</b>	10 x 10 cm	Box of 100 (50 pouches of 2 swabs)



## COTTON BALL

SKU	Description	Qty
<b>51630B</b>	Hydrophilic cotton balls, 0.7 g.	Bag of 700 pcs



## FIRST AID BANDAGE

- Dressing equipped with a non-adherent pad to the wound
- High absorption capacity
- Salmon-colored woven elastic support



SKU	Description	Dimension	Qty
EDM135102	Single dressing	7 x 2 cm	100

## EUROPLAIE BANDAGE DRESSING

- Sterile non-woven dressing
- Adhesive on all 4 sides with central absorbent gauze pad, covered with a Mylar film
- Does not adhere to the wound



- Suitable for light wounds or post-operative wounds with minimal exudate.

SKU	Description	Dimension	Qty
135261Z	Sterile dressing set	5 x 7,2 cm	100
135263Z	Sterile dressing set	10 X 10 cm	50

### Technical specifications:

- Class I sterile
- Ethylene oxide sterilization

## MICROPOROUS SURGICAL TAPE

- White non-woven backing
- High skin tolerance
- Pain-free removal



- To keep dressings on sensitive skin.
- Recommended for frequent dressing changes

SKU	Description	WD: with dispenser WoD: without dispenser	Dimensions	Qty
135300	Microporous surgical tape	WD	5 m x 2,5 cm	1 rolls
135325	Microporous surgical tape	WoD	5 m x 5 cm	6 rolls
135329	Microporous surgical tape	WD	5 m x 1,25 cm	24 rolls
135326	Microporous surgical tape	WD	5 m x 2,5 cm	12 rolls
135327	Microporous surgical tape	WD	9,14 m x 2,5 cm	12 rolls



## CLASSIC SUTURE SET

- Ready to use: saves time and efficiency
- Single-use set to prevent the risk of cross-contamination

### Set composition:

- 5 non-woven swabs 7.5 x 7.5 cm
- 1 sterile field 38 x 45 cm
- 1 anatomical forceps with teeth
- 1 pair of sharp scissors 11 cm
- 1 needle holder
- 1 rigid blister pack, 4 compartments

SKU	Description	Qty
EDMM50007/5	Suture kit	5



### Technical specifications:

- Class I sterile
- Ethylene oxide sterilization
- Single-unit packaging



## STERILE PLASTIC CUP

- For applying ointment

### Technical specifications:

- Class I sterile
- Ethylene oxide sterilization
- Single packaging
- Capacity: 120 ml
- Diameter: 80 mm

SKU	Description	Qty
EDM174617	Sterile plastic cup 120 ml	100



## MOLDED CELLULOSE SPITTOON

### À usage unique

- Disposable
- Disposable cellulose spittoon
- Biodegradable
- Water tightness of over 8 hours at a temperature of 40°C

### Technical specifications:

- Class I
- Capacity: 750 ml
- Weight: 22g

SKU	Capacity	Dimensions	Qty
125137A	Capacity 750 ml	253 x 146 x 44 mm	320



## ISOPROPYL ALCOHOL 70°

### Disinfectant and antiseptic solution

SKU	Description	Qty
201085	500 ml bottle	1
201107	250 ml bottle	1
208002	1L jerrycan	1

## CHLORHEXIDINE

### For fast and effective skin disinfection



SKU	Description	Qty
208003	50 ml spray bottle	1
208004	50 ml spray bottle	1

## SURGICAL DRAPE

### Single sterile

- Consists of an opaque polyethylene sheet lined with an absorbent cellulose sheet.

Scheme



- For minor surgery and special care

SKU	Description	Qty
161600	Single drape 40 x 45 cm	50
161601	Single drape 45 x 75 cm	50

### Technical specifications:

- Class I Sterile
- Sterilization with ethylene oxide

### Single sterile drape with adhesive

- Composed of an adhesive strip along its length, allowing the drape to be attached to another drape or to provide additional protection during the procedure.
- The drape is not likely to slip during the procedure.

Scheme




SKU	Description	Qty
161620	Sterile adhesive drapes 50 X 50 cm with 5 cm X 50 cm adhesive	50

### Fenestrated drape with sterile adhesive

- Combines the features of a fenestrated drape while ensuring its stability during the procedure.
- Diameter 6 cm + adhesive 12 x 12 cm around the hole.

Scheme

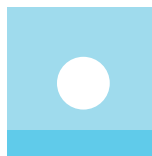

SKU	Description	Qty
161640	Sterile adhesive fenestrated drape 50 x 50 cm Diameter 6 cm + adhesive 12 X 12 cm	50
161641	Sterile adhesive fenestrated drape 50 x 70 cm Diameter 6 cm + adhesive 12 X 12 cm	50
161642	Sterile adhesive fenestrated drape 70 x 90 cm Diameter 6 cm + adhesive 12 X 12 cm	50

## SCALPELS

- Retractable safety scalpel
- Safety function before and after use to prevent accidental cuts
- Suitable for right and left-handed users
- Single-handed activation
- Audible position clicks
- Different sizes available in universal type 3 or 4 handle depending on blade size

### Handle type 3



### Handle type 4



#### Technical specifications:

- Class IIa
- Gamma ray sterilization
- Ribbed handle with engraved measuring scale from 0 to 5 cm
- Stainless steel blade
- ABS handle and guard

SKU	Description	Qty
<b>PNSF10</b>	Single-use sterile safety scalpel - size 10 - type 3 handle	10
<b>PNSF11</b>	Single-use sterile safety scalpel - size 11 - type 3 handle	10
<b>PNSF12</b>	Single-use sterile safety scalpel - size 12 - type 3 handle	10
<b>PNSF15</b>	Single-use sterile safety scalpel - size 15 - type 3 handle	10
<b>PNSF17</b>	Single-use sterile safety scalpel - size 17 - type 3 handle	10
<b>PNSF19</b>	Single-use sterile safety scalpel - size 19 - type 4 handle	10
<b>PNSF20</b>	Single-use sterile safety scalpel - size 20 - type 4 handle	10
<b>PNSF21</b>	Single-use sterile safety scalpel - size 21 - type 4 handle	10
<b>PNSF22</b>	Single-use sterile safety scalpel - size 22 - type 4 handle	10
<b>PNSF23</b>	Single-use sterile safety scalpel - size 23 - type 4 handle	10
<b>PNSF24</b>	Single-use sterile safety scalpel - size 24 - type 4 handle	10
<b>PNSF25</b>	Single-use sterile safety scalpel - size 25 - type 4 handle	10

## EUROFINE STANDARD BEVELED HYPODERMIC NEEDLE

- Various diameters and lengths available
- Very fine, offering smooth penetration without resistance, ideal for nearly painless injections
- Different sizes depending on the type of procedure and the patient's skin condition (thickness, fragility...)
- Great comfort for the patient



- To perform subcutaneous, intramuscular, intravenous injections, or samplings



SKU	Usage	Color	Gauge	Length	Diameter	Qty
101314	IM obese	Green	21 G	50 mm	0,8 mm	100
101306	IM Sampling	Green	21 G	40 mm	0,8 mm	100
101312	Vial puncture	Pink	18 G	40 mm	1,2 mm	100
101302	IV child	Blue	23 G	25 mm	0,6 mm	100
101301	Subcutaneous	Orange	25 G	16 mm	0,5 mm	100

### Technical specifications:

- Class IIa sterile
- Ethylene oxide sterilization
- Single use
- Sterile stainless steel needle
- Translucent polypropylene hub
- BN: Standard bevel
- Luer-type connection, compatible with all syringes on the market

## 3-PIECE SYRINGE

- 3-piece syringe: piston + elastomeric seal + body
- Double-contact sealing point
- Centered Luer Lock tip
- Graduated in ml

SKU	Description	Capacity	Qty
103330C	3-piece syringe	5 ml	100
103332c	3-piece syringe	10 ml	100
103334C	3-piece syringe	20 ml	100
103410	3-piece syringe	50 ml	25

### Technical specifications:

- Class IIa sterile
- Ethylene oxide sterilization



## SHARPSAFE WASTE CONTAINER

- Designed for collecting sharp medical waste such as scalpel blades, needles, syringes...
- Made from recycled plastic materials

SKU	Description	Capacity	Qty
162504	Waste container	4L	50
162505	Waste container	7L	50

### Technical specifications:

- Sharpsafe®, certified EN ISO 23907-1:2019 (mandatory in France since 01/01/2022)
- NF Mark validation





## NUCLEAN®

### Spray & wipes

ALCOHOL-FREE !



#### Total activity in 5 minutes

- **Two versions for 2 uses:**

NuClean Spray 500 ml & NuClean wipes: Disinfection by wiping of equipment (monitors, screens) and medical devices (stretchers and examination tables).

- **NuClean SH Spray 750 ml:**

Disinfection by wiping of all high surfaces in medical environment (countertops, carts, and furniture). Alcohol-free.

**Microbiological properties:**

	Spray & lingettes NuClean	Spray NuClean SH
<b>BACTERICIDE</b>	EN 13727: active in 1 min	EN 13697 : active in 1 min. EN 13727 : active in 5 min. NFT 72-281 : active in 5 min.
<b>YEASTS/ MOLDS</b>	EN 16615 (dirty conditions with mechanical action): Pseudomonas aeruginosa, Staphylococcus aureus, Enterococcus hirae: active in 5 min	EN 13624 , EN 1650, EN 13697, NF T 72-281: active in 5 min.
<b>VIRUCIDE</b>	EN 14476: active in 5 min.	EN 14476 , EN 16777: active in 5 min.
<b>SPORICIDE</b>	EN 13704: active in 15 min.	NF T 72-281: active in 5 min.

Find the full microbiological properties on the technical data sheet available upon request.

SKU	Description	Qty
AL4004	Dispenser box of 120 wipes - 130 x 190 mm	1 unit
EC500	500 ml spray bottle	1 unit
EC750SH	750 ml spray bottle - Alcohol-free for high surfaces	1 unit

## SURFA' SAFE

### Spray

- For cleaning and disinfecting surfaces
- Fast antimicrobial efficacy: bacteria including MRSA, Rotavirus... in just 2 minutes.



**Microbiological properties:**

**BACTERICIDE:** EN 13727, EN 13697: active in 2 min  
**YEAST / FUNGICIDE:** EN 13697: active in 5 min, EN 13624: active in 2 min  
**MOULDS:** EN 1275, EN 13624, EN 13697: active in 20 min  
**VIRUCIDE:** EN 14476: VIH, PRV, BVDV, Vaccine, HSV, Rotavirus: 1 min  
 EN 14476: Polyomavirus SV 40, Coronavirus: 5 min  
 EN 14476: Adenovirus, Norovirus murin: 30 min

Find the full microbiological properties on the technical data sheet available upon request.

SKU	Description	Qty
AN842326	750 ml spray bottle	1 unit

## ANIOS MILD SOAP

### High frequency

- Simple hand washing, specially formulated for sensitive skin and repeated use
- Contains glycerin

SKU	Description	Qty
AN205228	500 ml bottle with pump	1 unit
AN205229	1L bottle with pump	1 unit
AN205034	5L jerrycan	1 unit



## ANIOSGEL 800

- Hydro-alcoholic gel for hygienic treatment and surgical disinfection (EN12791) of hands
- Can be used without water or rinsing.
- Fragrance-free, hypoallergenic formula

### Microbiological properties:

**BACTERICIDE:** EN 1040, EN 13727  
(dirty conditions: albumin and erythrocytes MRSA (EN 13727): 5 minutes.  
**MYCOBACTERICIDE:** EN 14348 : 30 sec  
**LEVURICIDE:** EN 1275, EN 13624 (2013): 30 sec  
**VIRUCIDE:** EN 14476 (2013) on Poliovirus : 60 sec ; Adenovirus et Norovirus Murin: 30 sec. According to EN 14476 (2013) on BVDV (HCV model virus), PRV (HBV model virus), HIV-1, HSV-1, Rotavirus, RSV, Vaccinia virus: 30 sec.

SKU	Description	Qty
AN2550505	100 ml bottle with pourer spout	1 unit
AN2550762	300 ml bottle with pump	1 unit
AN2550748	500 ml bottle with pump	1 unit
AN1644679	ANIOSGEL85 NPC Refill - 1L clear, fragrance-free for healthcare settings	1 unit
AN425163	Anios wall-mounted automatic dispenser for healthcare settings	1 unit



## HAND TOWELS & TOILET PAPER

SKU	Description	Qty
708	Interfolded white 2 plies tissue napkins (packed in 30 packs of 125), unfolded 22 x 35 cm, folded 22 x 9 cm	1 unit
709	Hand towels 25 x 38 cm, lateral dispensing	24
710	Hand towels 22 x 38 cm, central dispensing	12
S401023	Toilet paper, 200 sheets/roll	96



# EQUIPMENT

- Medical furniture
- Office supplies
- Diagnosis



## EXAMINATION CHAIR



### LE QUESNOY

### 2006/1

<b>Type of chair</b>	3 electric sections	3 sections
<b>Adjustable parts</b>	Leg rest and backrest	Backrest
<b>Backrest / leg rest control</b>	<ul style="list-style-type: none"> <li>• Electric</li> <li>• Gas spring</li> </ul>	By gas spring
<b>Possible positions</b>	<ul style="list-style-type: none"> <li>• Gynecological examination</li> <li>• Ultrasound (Horizontal position)</li> <li>• Trendelenburg position (Trendelenburg/ reverse Trendelenburg)</li> </ul>	<ul style="list-style-type: none"> <li>• Gynecological examination</li> <li>• Ultrasound (Horizontal position)</li> <li>• Trendelenburg position (Trendelenburg/reverse Trendelenburg) - Optional</li> </ul>
<b>Variable height</b>	Variable height adjustment from 620 to 1020 mm	Fixed
<b>Type of manual control</b>	Yes	-
<b>Adjustable headrest</b>	Yes	No
<b>Gynecological stirrups</b>	Yes	Yes
<b>Pair of fixed foam stirrups</b>	Included as standard	Optional
<b>Pair of stainless steel stirrups</b>	Included as standard	Optional
<b>Removable tray under the seat</b>	Yes	Yes
<b>Roll holder on the backrest</b>	Yes	Yes
<b>Dimensions in cm (L x W x H)</b>	175 x 80 x 64 à 104	162 x 70 x 80
<b>Chair weight</b>	90 kg	45 kg
<b>Maximum weight</b>	300 kg	220 kg
<b>Seat width in</b>	60	60
<b>Manufacturer's warranty</b>	3 years	3 years
<b>Type of braking</b>	4 casters with Ø 100 mm brakes or with centralized wheel locking.	Fixed chair
<b>Options</b>	<ul style="list-style-type: none"> <li>• Pair of polyurethane stirrups with ball joint system. (<b>ROTULES/PU</b>)</li> <li>• Transparent PVC crystal protection (<b>PROTEC</b>)</li> <li>• IV pole support with 2 single hooks (<b>3511/0</b>)</li> <li>• Notch on the seat of the couch (<b>ECHANC/3PLANS</b>)</li> </ul>	<ul style="list-style-type: none"> <li>• Pair of stainless steel stirrups. (<b>1991/0</b>)</li> <li>• Pair of polyurethane foam stirrups. (<b>1991/2</b>)</li> <li>• Pair of PU stirrups with ball joint system. (<b>ROTULES/PU</b>)</li> <li>• Assise relevable (position de Trendelenburg -5°/+12°). (<b>2006/2</b>)</li> <li>• Retractable footrest (<b>3096/1</b>)</li> <li>• Notch on the seat of the couch. (<b>ECHANC/3PLANS</b>)</li> </ul>

## STOOL

Available colors (only for diameter 360 mm):



SKU	Description
3600	Stool
3700	Backrest option

### Technical specifications

- Seat color: diameter 360 mm, thickness 60 mm
- Upholstery classified M1 or M2
- Chromed gas spring
- Adjustable height from 530 mm to 720 mm
- 5-star base
- 50 mm diameter casters



## STEP STOOL

SKU	Description
2010	Step stool 2 steps
2000	Step stool 1 step

### Technical specifications

- Maximum weight 110 kg
- Frame made of stainless steel tube
- Welded structure

2000



2010



## TROLLEY

SKU	Description	Dimensions
55	Stainless steel two-tier trolley without handle, color	600 x 400 mm

### Technical specifications

- Trolley fully welded stainless steel structure
- Frame with a shiny polished 25 mm diameter
- Composite resin color trays
- 100 mm diameter ABS wheels with bumper rollers





Medical furniture



EDM220



EDM320

Technical specifications

- Trolley structure in color composite resin and stainless steel tubes with a diameter of 30 mm
- Smooth composite resin top trays without guiding handles
- 100 mm diameter ABS wheels with bumper rollers
- Available in 2 or 3-tier versions

Options

- Drawers
- Dividers
- Rails
- Waste bag holders

SKU	Description	Dimensions
220	Compact two-tier trolley	650 x 450 mm
320	Compact three-tier trolley	650 x 450 mm

Please specify your color when placing your order



SERVICE TROLLEY



Technical specifications

- Fully welded stainless steel structure without screws
- Width of the 3 panels: 54, 42, and 54 cm
- 50 mm diameter polyamide wheels
- PVC-coated polyester textile panels
- Fire resistance classification M1
- Removable and washable canvas

SKU	Description	Dimensions
312	Service trolley	790 x 480 x 950 mm

SCREEN



Technical specifications

- Fully welded stainless steel structure without screws
- Width of the 3 panels: 54, 42, and 54 cm
- 50 mm diameter polyamide wheels
- PVC-coated polyester textile panels
- Fire resistance classification M1
- Removable and washable canvas

Please specify your color when placing your order



SKU	Description	Dimensions W x H
8330	Stainless steel mobile screen	129 x 171 cm

## MEDICINE CABINET

### Technical specifications

- Body made of 8/10mm white sheet metal. Swing door(s) with a 2-point lock.
- Supplied with 4 glass shelves.
- 1 Door - weight: 54 kg.
- 2 Doors - weight: 75 kg.

SKU	Description	Dimensions
<b>SML101</b>	Single-door medicine cabinet	600 x 1800 x 400 mm
<b>SML102</b>	Double-door medicine cabinet	900 x 1800 x 400 mm



## WALL CABINET

### Technical specifications

- Body made of 8/10mm white sheet metal only.
- Equipped with 1 door with magnetic closure, double magnet, and 2 glass shelves.
- Weight: 10 kg.

SKU	Description	Dimensions (W X H X D)
<b>SZL101</b>	Single-door Wall cabinet	400 x 500 x 250 mm



### Technical specifications

- Body made of 8/10mm white sheet metal only.
- Equipped with 2 sliding glass doors and 2 glass shelves.
- Weight: 20 kg.

SKU	Description	Dimensions (W X H X D)
<b>SZL102</b>	Double-door Wall cabinet	1000 x 500 x 250 mm



## GYNECOLOGY LAMP

- LED technology without infrared: reduces heat sensation for both patient and doctor and saves energy
- Positioning and long r1 unit with flexible arm that avoids blind spots
- Excellent stability
- Dimmer located on the control box
- Ingeniously placed switch on the lamp head

### Technical specifications

- Illuminance: D = 30 cm 80,000 Lux, D = 50 cm 35,000 Lux
- Color temperature: 5000-6000K
- Average LED lifespan: over 35,000 hours
- Power consumption: 85-265V 50/60Hz
- Type of extension: 76 cm gooseneck
- Rod height: 71 cm

SKU	Description
<b>OSLAMP01</b>	Black examination lamp
<b>OSLAMP02</b>	White examination lamp



## BATTERY



### Rechargeable batteries

SKU	Description	Qty
635177	Energizer Rechargeable Battery LR3 (AAA) NiMH Power Plus HR03 700 mAh 1.2V	4
635730	Energizer Recharge Extreme AA 2300 mAh	4



### Lithium batteries

SKU	Description	Qty
CR2025L/1BP	Panasonic CR2025 Lithium Battery	10
CR-2032EP/IB	Panasonic CR2032 3V Lithium Button Cell Battery	12
628287	Energizer 2CR5 6V Lithium Battery	1
635233	Ultimate L92 AAA 1.5V Lithium Battery	48



### Piles Alcaline

SKU	Description	Qty
LR03AD/4P	Panasonic LR03 AAA 1.5V alkaline battery	48
LR14AD/AP	Panasonic LR14 1.5V alkaline battery	24
LR20AD/4P	Panasonic LR20 1.5V Alkaline battery	24
LR44L/1BP	Panasonic LR44 alkaline battery	12
LR6AD/4P	Panasonic LR06 1.5V alkaline battery	48
3LR12PPG/1BP	Pro power 4.5V alkaline battery	12
6LR61AD/1P	Panasonic 6LR61 9v alkaline battery	12
608305	Energizer 3LR50 alkaline battery	10

# Supporting you every day

BIOPSY



IMAGING



DISINFECTION



GYNECOLOGY



MEDICAL FURNITURE



RADIATION PROTECTION



PROTECTIVE COVERS



PROBE COVER



NON-MAGNETIC FURNITURE



(REF: C-FEMINA-EN-EXP. CE04-2024)

*femina*  
by EDM

**EDM**  
medical  
imaging

**EDM**

1, rue des Boisseliers - 95330 Domont - France  
[www.edm-imaging.com](http://www.edm-imaging.com)

Tél.: + 33 (0)1 30 11 83 50

Fax: + 33 (0)1 39 35 75 25